

AD-1688 Specifications

Applicable devices	Balances and scales with an RS-232C connector manufactured by A&D
Memory size	Approx. 5,000 data sets in the standard format (15 characters in 1 data set) can be stored
Clock error	±1 minute / month
Applicable OS	Windows 2000 / XP / Vista
Dedicated power supply	Not required. Power is supplied by the weighing instrument or PC.
Operating environment	5 °C to 40 °C (41 °F to 104 °F), 85 %RH or less (no condensation)
Dimensions	55 × 103 × 16.5 mm
Weight	Approx. 60 g

Standard Accessories

USB extension cable (0.5 m)
Three RS-232C cables (1 m)
 with a D-sub 9-pin connector
 with a D-sub 25-pin connector
 with a 7-pin DIN connector

The data is stored in CSV format and can be used easily by software such as Microsoft Excel.

	A	B	C	D	E	F
1	1	6/10/2009	10:41:24	ST	3.076	g
2	2	6/10/2009	10:41:36	ST	3.082	g
3	3	6/10/2009	10:42:06	ST	3.102	g
4	4	6/10/2009	10:42:38	ST	3.058	g
5	5	6/10/2009	10:42:55	ST	3.091	g
6	6	6/10/2009	10:43:01	ST	3.089	g

Quick USB Adapter AD-8527




Transmits weighing data directly from an A&D balance or scale to a USB port of a PC (no driver software is required)*

Sends weighing data to software applications such as Microsoft Excel, Word, and Notepad

Small size and IP65 dust and waterproof (with the casing cover)

*Only unidirectional communication from a balance or scale to a PC is possible.

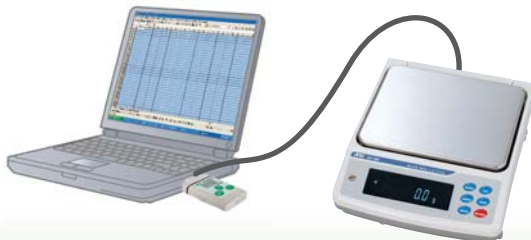

AD-8527 Specifications

Applicable devices	Balances and scales with an RS-232C connector manufactured by A&D
Dedicated power supply	Not required. Power is supplied by the weighing instrument or PC.
Operating environment	5 °C to 40 °C (41 °F to 104 °F), 85% RH or less (no condensation)
Dimensions	55 × 103 × 16.5 mm
Weight	Approx. 60 g
Standard accessories	USB extension cable (0.5 m), Three RS-232C cables (1 m)

A&D

...Clearly a Better Value

A&D Company, Limited

3-23-14 Higashi-Ikebukuro, Toshima-ku, Tokyo 170-0013 JAPAN
 Telephone:[81](3) 5391-6132 Fax:[81](3) 5391-6148
<http://www.aandd.jp>

A&D ENGINEERING, INC.

1756 Automation Parkway, San Jose, CA 95131 U.S.A.
 Telephone:[1](408) 263-5333 Fax:[1](408) 263-0119

A&D Australasia Pty Ltd.

32 Dew Street, Thebarton, South Australia 5031 AUSTRALIA
 Telephone:[61](8) 8301-8100 Fax:[61](8) 8352-7409

A&D INSTRUMENTS LTD.

Unit 24/26 Blacklands Way Abingdon Business Park,
 Abingdon, Oxon OX14 1DY UNITED KINGDOM
 Telephone:[44](1235) 550420 Fax:[44](1235) 550485

<German Sales Office>

Große Straße 13 b 22926 Ahrensburg GERMANY
 Telephone:[49](0) 4102 459230 Fax:[49](0) 4102 459231

A&D KOREA Limited

Manhattan Bldg. 8F, 36-2 Yoido-dong, Youngdeungpo-gu, Seoul, KOREA
 Telephone:[82](2) 780-4101 Fax:[82](2) 782-4280

A&D RUS CO., LTD.

Vereyskaya str.17, Moscow, 121357 RUSSIA
 Telephone: [7] (495) 937-33-44 Fax: [7] (495) 937-55-66

A&D Instruments India Private Limited

509 Udyog Vihar Phase V
 Gurgaon-122 016, Haryana, INDIA
 Telephone: [91](124) 471-5555 Fax: [91](124) 471-5599